

Roho[®] Enhancer[®] Cushion

Individual air cells allow
air circulation between
the cells

Durable Neoprene, Flame
resistant and resistant to
urine, oils and medicines

No user
weight limit



Two-manifold system for enhanced
midline channeling of the femurs,
lateral stability and tissue protection

Stability

The Roho[®] Enhancer[®] Cushion provides the skin protection you expect from a Roho Group product, with the positioning and stability needed for increased function.

Positioning

By offering a contoured form with a two manifold system, the client/user can benefit from positive midline channelling of the legs along with increased lateral stability.

Skin Protection

In addition, the special design of the Enhancer cushion inhibits sacral sitting while enhancing tissue viability.



Head Office: 2/19 Boden Road, Seven Hills, NSW 2147

Phone: 02 9674 8904

Fax: 02 8212 5912

Website: www.dmedirect.com.au

Roho® Enhancer® Cushion

The Roho® Enhancer® Cushion offers you a unique option in wheelchair seating. By offering a contoured form with a two manifold system, the client/user can benefit from positive midline channelling of the legs along with increased lateral stability. In addition, the special design of the Enhancer cushion inhibits sacral sitting while enhancing tissue viability.

The Enhancer unique contoured shape is the result of flexible air cells of varying heights. The two manifold system allows you to: 1) adjust the interior bucket seat section for reduced pressure on the ischia and coccyx while 2) adjusting the two outer perimeter rows for enhanced stability and positioning.

The two design features combined offer you a cushion that shapes to the user's body providing superior pressure distribution and minimal shear forces that are characteristic of every Roho® Dry Floatation® cushion, while offering you the added benefit of progressive position.

FEATURES AND BENEFITS

- Adjustable:** Fits , matches and tracks client shape to manage integrity, positioning and stability
- Lightweight:** Adds very little weight to mobility system - 1.8kg
- Washable:** Can be easily cleaned with soap and water. Cover is machine washable.
- Repairable:** Neoprene rubber cushion can be patched or repaired to extend use.
- Ventilated:** Separate air cells allow air circulation to reach surface of the skin and provide a cooling environment.

SPECIFICATIONS

Includes: Enhancer Cushion, 2 way stretch cover, hand inflation pump, repair kit, operating instructions

Construction: Flame resistant, neoprene rubber


Sizes: Range from 37.5cm up to 50cm

User Weight Limit: NONE

MORE PRODUCTS



Jet Stream Pro® Back Support System
Ultra lightweight back rest

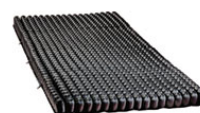

ROHO® Dry Floatation® Mattress Overlay System