


# Roho<sup>®</sup> Hybrid Elite<sup>®</sup> Cushion

## Dual Compartment Cushion

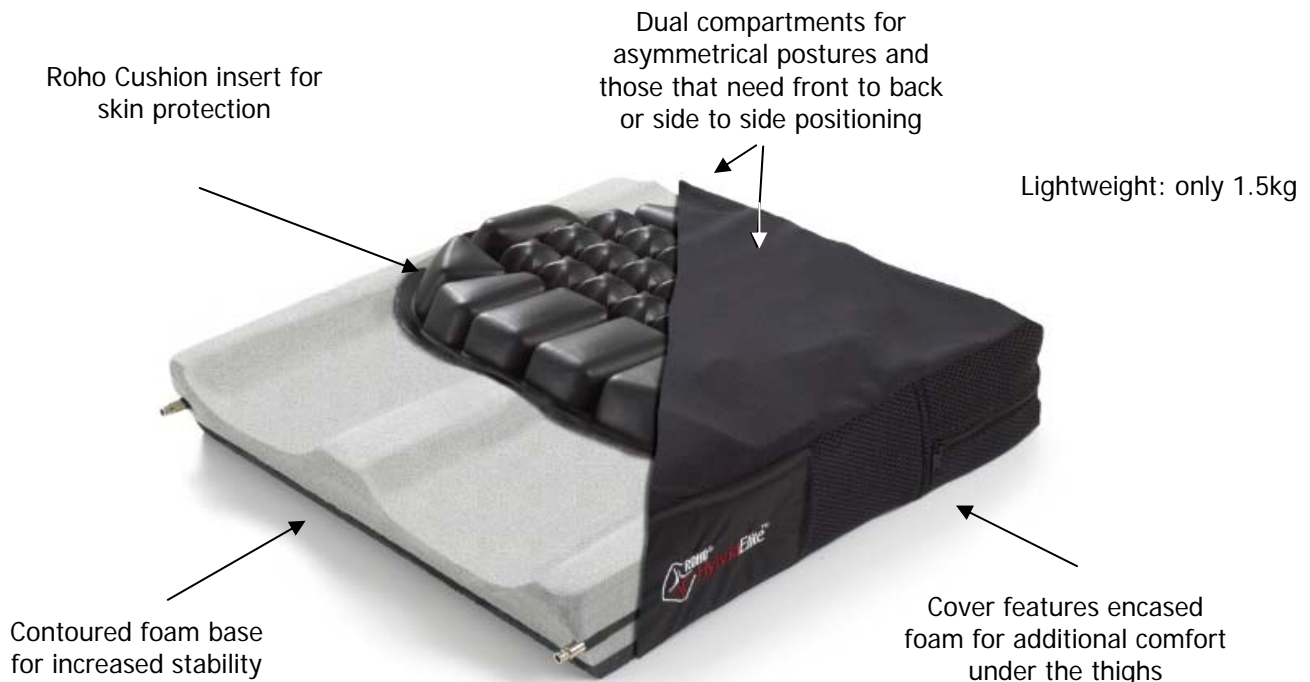

### Positioning

### Skin Protection

### Lightweight

The Roho<sup>®</sup> Hybrid Elite<sup>®</sup> Cushion Dual Compartment Cushion combines the superb skin protection and shape matching ability of Roho Dry Floatation<sup>®</sup> with the stability of Jay<sup>®</sup> base technology.

New cushion design allows for improved pelvic alignment and stability for improved position and function while maintaining optimum immersion and envelopment of both the ischial tuberosities and greater trochanters.

Available up to 60cm seat width with a 225kg user weight limit, ideal for Bariatric Clients.

# Roho® Hybrid Elite® Cushion Dual Compartment Cushion

## Positioning


Achieve greater ability to position your client with the Roho® Hybrid Elite® Cushion Dual Compartment Cushion. Its customized Jay contoured foam base provides you with the opportunity for both lower extremity and pelvic alignment and control.

## Skin Protection

The Hybrid Elite Dual Compartment Cushion hydrostatically loads the ischial tuberosities and the posterior shelf of the greater trochanters for those clients who are at high risk of breakdown. Unparalleled in performance, the ROHO cushion insert allows for immersion and envelopment of the client's individual shape. It provides an adjustable, reduced friction and shear environment that can help protect the client from acquiring pressure ulcers from sitting. The dual valve option on the ROHO insert allows for additional customization to match a client's specific shape and asymmetry.

## Stability

The Hybrid Elite Dual Compartment Cushion's customized Jay contoured foam base provides improved positioning of the pelvis and lower extremities, enhancing sitting posture. This added stability assists the client in achieving a

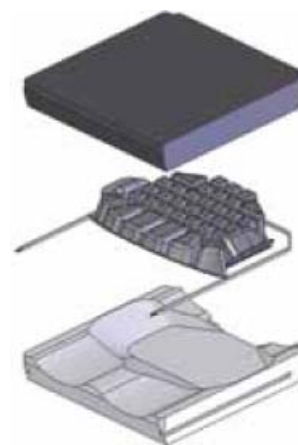

Roho Hybrid Elite Dual Compartment Cushion HE-2RHYxxyyC

## FEATURES AND BENEFITS

- Adjustable:** Cushion insert matches the individual's shape by adjusting the amount of air.
- Washable:** Air cushion insert is easily cleaned with soap and water. Cover with encased comfort foam is machine washable. Foam base may be wiped with damp cloth.
- Repairable:** Neoprene rubber cushion insert can be patched or repaired to extend use.
- Fluid Resistant:** Customized Jay foam base and ROHO air insert are fluid resistant.

## SPECIFICATIONS

### Construction:

- ROHO cushion insert:  
flame resistant, neoprene rubber
- Customized Jay contoured foam base:  
flame resistant, fluid resistant closed cell foam

**Sizes:** Range from 35cm up to 60cm

**User Weight Limit:** 225kg

## MORE PRODUCTS



**Elevation**  
Ultra lightweight back rest



**Jet Stream Pro® Back Support System**  
Ultra lightweight back rest