


Roho[®] Quadtro Select[®] High Profile[®] Cushion

5cm cell height offers therapeutic performance for the treatment and prevention of ischemic ulcers

Individual air cells allow air circulation between the cells

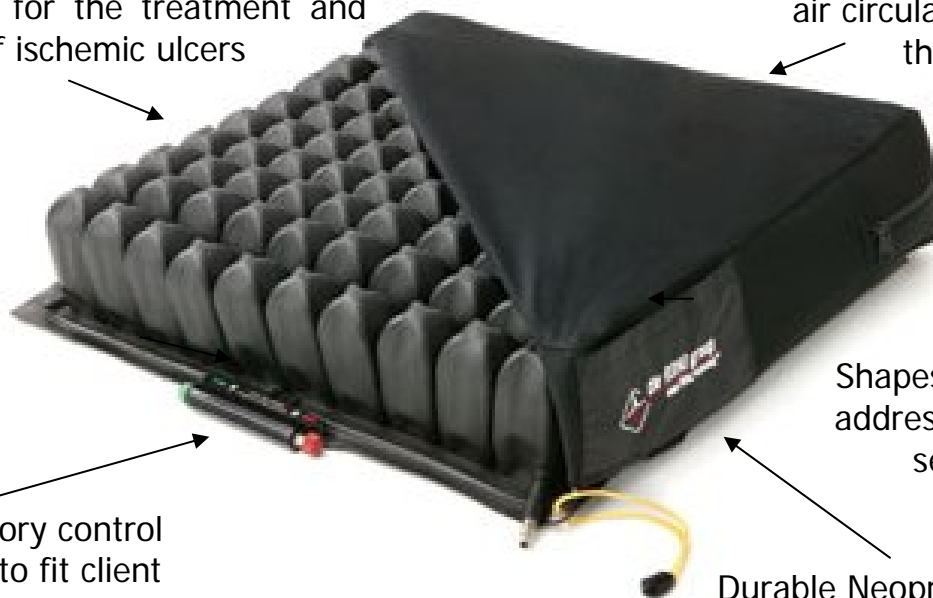

Isoflo memory control allows you to fit client while seated

No user weight limit

Shapes to the individual addressing their specific seating needs

Durable Neoprene, Flame resistant and resistant to urine, oils and medicines

Positioning

Skin Protection

Stability

Convenience

The Roho[®] Quadtro Select[®] High Profile[®] Cushion sets the standard in the industry for overall performance in wheelchair seating.

Each and every user will benefit from a customised fit through the simple push of a knob. The revolutionary ISOFLO Memory Control offers shape fitting capabilities while the user is seated, allowing quick and easy, on-demand adjustment to maximise function.

With the built in stability and simplicity of the Quadtro Select High Profile Cushion, no longer will you have to sacrifice maximum skin protection to get stability, positioning or convenience. The Quadtro Select High Profile Cushion does it all. And with no compromise.


Roho[®] Quadtro Select[®] High Profile[®] Cushion

The Roho[®] Quadtro Select[®] High Profile[®] Cushion allows for quick, simple position selection for every user while seated. By pushing the knob to the open position, you can adjust the users posture to the desired sitting position. By pushing the knob to the closed position you can maintain the desired sitting posture for maximum function.

Based on thirty years research and client experience patented Roho[®] Dry Floatation[®] technology provides unparalleled performance in facilitating blood flow by fitting matching and tracking the shape of a person.

The Quadtro Select High Profile Cushion offers a new level of stability never before found in an air floatation cushion. Stability is enhanced by locking air into place in each of the four quadrants through the simple use of the isoflo memory control, minimising side to side or front to back motion.

All postural adjustments are made while the user is seated, eliminating the need for transfer on and off the chair during the evaluation process. In addition the Quadtro Select High Profile Cushion allows the user to selectively change their position throughout the day as their particular activity demands.

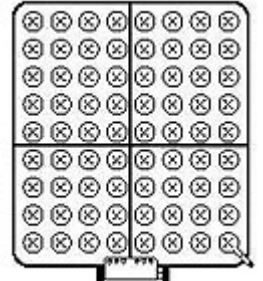

FEATURES AND BENEFITS

- Adjustable:** Fits specific seating needs and unique body shapes
- Lightweight:** Easy to handle and carry at just 1.7kgs
- Washable:** Can be easily cleaned with soap and water. Cover is machine washable
- Repairable:** Neoprene rubber cushion can be patched or repaired to extend use
- Ventilated:** Separate air cells allow air circulation to reach surface of the skin and provide a cooling environment

SPECIFICATIONS

Includes: Quadtro Select Cushion, 2 way stretch cover, hand inflation pump, repair kit, operating instructions

Construction: Flame resistant, neoprene rubber


Sizes: Range from 30cm up to 55cm

User Weight Limit: NONE

MORE PRODUCTS



Jet Stream Pro[®] Back Support System
Ultra lightweight back rest

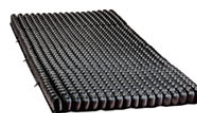

ROHO[®] Dry Floatation[®] Mattress Overlay System