


Hand Splints

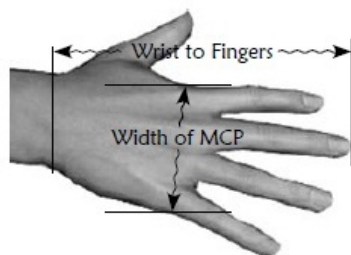

| SIZE | WIDTH OF MCP | WRIST TO FINGER |
|---------|--------------|-----------------|
| X-Small | 7cm - 7.5cm | Up to 13cm |
| Small | Up to 7.5cm | Up to 16.5cm |
| Medium | Up to 9cm | Up to 18cm |
| Large | Up to 9.5cm | Up to 20cm |

HYPER HAND PAEDIATRIC

MR20717



- Uses:** Treats from the tightest closed palm, to moderate flexion.
- Features:** Soft silky padding and straps. Plastic slit between each digit can be moulded per finger.
- Sizes:** PEDIATRIC, Infant, Toddler, Child, Youth

RESTORATIVE HAND

MR20736



- Uses:** Effectively and comfortably treat any wrist/hand/finger lost range and tension from fresh stroke, brain or spinal cord injuries to moderate deformities and pain.
- Features:** Flex design encourages joints to relax. Heat mould to provide just low load passive stretch on the joints, then remould to gradually correct joints to normal alignment. Provides prolonged low load passive stretch to re-lengthen shortened tissue. Kydex base extends to end of thumb. Comfortably supports Osteo and Rheumatoid arthritic joints and treats. Prevents deformities if treated in early stages.
- Sizes:** ADULT Small, Medium, Large
- Options:** 2 attachments - the Dorsal Interphalangeal Assist and the Dorsal Attachment

NEUROFLEX FLEX HAND

MR20720 Adult
MR20720P Ped



- Uses:** Treat from the tightest closed palm, to moderate flexion.
- Features:** Efficiently works to relax joint tension, even from the wrist. Graduated sizes of cones gradually open hand/fingers/thumb.
- Sizes:** PEDIATRIC Universal (One Size); ADULT Universal (One Size)

NEUROFLEX FLEX THUMB EASE


MR20737P Ped
MR20737 Adult



- Uses:** Ideal for severe thumb deformities.
- Features:** Gradually add more filler to pocket under thumb to slowly bring it to better alignment. Mould Kydex base to apply only low load passive stretch, then remould to gradually take fingers/wrist from deformity toward correction. Comfortable. Completely covered with washable, breathable, moisture wicking padding with removable finger separators.
- Sizes:** PEDIATRIC Infant, Toddler, Child, Youth ADULT Universal (One Size)

Please specify Left or Right when ordering

Hand Splints

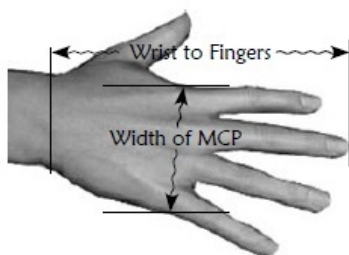

| SIZE | WIDTH OF MCP | WRIST TO FINGER |
|---------|--------------|-----------------|
| X-Small | 7cm - 7.5cm | Up to 13cm |
| Small | Up to 7.5cm | Up to 16.5cm |
| Medium | Up to 9cm | Up to 18cm |
| Large | Up to 9.5cm | Up to 20cm |

DORSAL INTERPHALANGEAL FLEXION ASSIST

MR20736D



- Uses:** For extensor tone or hyper extended joints.
- Features:** Elastic straps placed over ends of fingers, provides non-rigid treatment.

DORSAL ATTACHMENT

MR20736I



- Uses:** Addresses extensor tone in the hand.
- Features:** Semi-rigid component turns any of our hand splints into a dorsal splint. Place it on the back of the hand secured with the hand splint's straps. It is an additional splint along top of arm. Remould the dorsal attachment and the base hand splint simultaneously to "sandwich" the hand in between.
- Sizes:** ADULT & Pediatric

THUMB SPICA

MR20719



- Uses:** Treats injuries like strains, sprains, & repetitive motion injuries like Dequeveins Syndrome. Excellent for patient compliance.
- Features:** Immobilizes thumb & wrist comfortably while leaving fingers free. Supports while joints heal. Remouldable plastic is totally covered with moisture wicking, washable padding.
- Sizes:** ADULT Small, Medium, Large

COCK-UP SPLINT

MR20739



- Uses:** Treat Carpal Tunnel Syndrome & other conditions that require stabilizing the wrist.
- Features:** Holds wrist in desired degree of cock-up position/shipped at 30 degrees. Easily remoulded with heat. Comfortable with soft strap through web space. Washable lining.
- Sizes:** ADULT Small, Medium, Large

Please specify Left or Right when ordering