


Spinal Braces

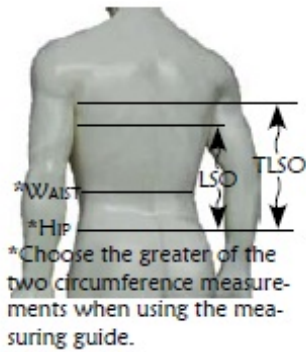

PEDIATRIC

SIZE	WAIST/HIP CIRC	TLSO HEIGHT	FRONT HEIGHT
Infant	46cm - 63.5cm	18cm	10cm
Toddler	46cm - 63.5cm	25.5cm	15cm
X Short	63.5cm - 76cm	30.5cm	18cm

ADULT

SIZE	WAIST/HIP CIRC	SHORT LSO	AVE. LSO	TALL LSO
XSLIM	63.5cm - 76cm	30.5cm	33cm	35.5cm
SLIM	76cm - 89cm	30.5cm	33cm	35.5cm
REGULAR	89cm - 101.5cm	30.5cm	33cm	35.5cm
FULL	101.5cm - 117cm	30.5cm	33cm	35.5cm
XFULL	117cm - 134.5cm	30.5cm	33cm	35.5cm

VERTEBREASE TLSO

MR20730



- Uses:** Treat neurological tone, lost range of motion, chronic low back pain, spondylolysis, degenerative disc disease, facet syndrome, post-op fusion, spondylolisthesis, compression fracture, spinal fractures, and spinal stenosis.
- Features:** Four overlapping Kydex® plastic shells restricts gross trunk movement for tri-planar control. Thermoplastic allows for customized heat moulding. Shoulder straps improve posture and relieve pain.
- Sizes:** PAEDIATRIC (See Sizing Chart); ADULT (See Sizing Chart)

VERTEBREASE LSO

MR20729



- Uses:** Treats conditions like compression fractures, osteoporosis, degenerative disc disease, spondylolisthesis, spinal stenosis, and also Kyphosis (upper body forward leaning).
- Features:** Appropriate for Orthopedic or Restorative patients. Effectively treat obese or pregnant patients or those with uniquely shaped abdomens with the optional pendulous abdomen front. Shoulder straps improve posture and relieve pain. Cervical attachments available to correct head/neck alignment - Hyperextension Stop, Safe Cerv Left or Right and Cervical Extension
- Sizes:** ADULT (See Sizing Chart)

SAFE SPINE TLSO

MR20748



- Uses:** Comfortably treat functional curves of spine/back like Scoliosis, Lordosis, and Kyphosis. Supporting the spine and relieving pain from common conditions like osteoporosis and degenerative disc disease.
- Features:** Prevent deformities of the spine by supporting patients with decreased or lack of muscle tone as in Cerebral Palsy and certain brain and spinal cord injuries. Shoulder straps encourage shoulders and upper body toward and to corrected posture. Helps to relieve pain. Wear while walking or sitting. Cervical attachments available to correct head/neck alignment - Hyperextension Stop, Safe Cerv Left or Right and Cervical Extension.
- Sizes:** ADULT (See Sizing Chart)

Spinal Braces

CERVICAL EXTENSION

MR20728C



- Uses:** Treats forward head/upper body posture.
- Features:** Attaches to back of VertebrEase LSO, TLSO or Safe Spine TLSO. Comfortable, non rigid. Heat mouldable for best fit, and remould to gradually improve alignment. Use while walking or sitting.
- Sizes:**

HYPER-EXTENSION

MR20728H



- Uses:** Treats hyper extended head leaning.
- Features:** Attaches to back of VertebrEase LSO, TLSO or Safe Spine TLSO. Heat mouldable for exact fit. Remould/reset to gradually correct alignment. Semi-rigid design is comfortable & applies prolonged low load passive stretch to reach relaxed state and re-lengthen shortened tissue, correcting deformities.
- Sizes:**

SAFECERV L/R

MR20728C-L/R



- Uses:** Treats left head leaning with or without forward leaning.
- Features:** Attaches to back of VertebrEase LSO, TLSO or Safe Spine TLSO. Heat mouldable for exact fit. Remould/reset to gradually correct alignment. Semi-rigid design is comfortable & applies prolonged low load passive stretch to reach relaxed state and re-lengthen shortened tissue, correcting deformities.
- Sizes:**